
waterpik®

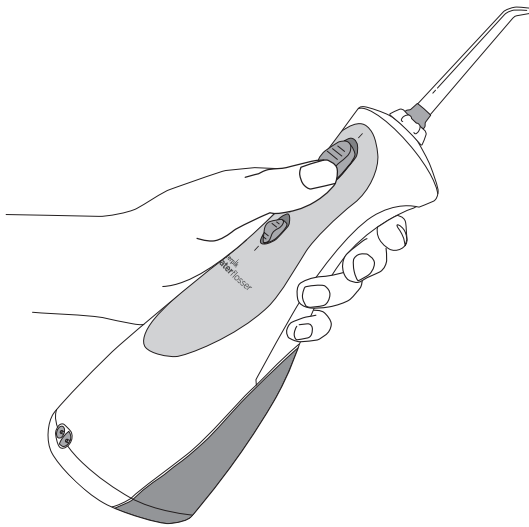**Fio Dental de Água
Waterpik®
Modelo WP-450/462**

**Irrigador Oral
Waterpik®
Modelo WP-450/462**

**Waterpik®
Water Flosser
Model WP-450/462**

.....
www.waterpik.com

samcooper**design**

DO NOT PRINT

DATE: 04.26.16

CLIENT: WATERPIK

FILENAME: 20021794-FAA_WP450B_IM_revs_rd6.indd

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Important Safeguards	21
Product Description	23
Getting Started	25
Cleaning and Troubleshooting Your Waterpik® Water Flosser	26
Limited One-Year Warranty	27

IMPORTANT SAFEGUARDS

When using electrical products, especially when children are present, basic safety precautions should always be followed, including the following:

READ ALL INSTRUCTIONS BEFORE USING.

DANGER:

To reduce the risk of electrocution:

- Do not handle charger with wet hands.
- Do not place in or drop into water or other liquid.
- Do not use while bathing.
- Do not place or store product where it can fall or be pulled into a tub or sink.
- Do not reach for a product that has fallen into water. Unplug immediately.
- Check the charger cord for damage before the first use and during the life of the product.


WARNING:

To reduce the risk of burns, electrocution, fire, or injury to persons:


- Do not plug this device into a voltage system that is different from the voltage system specified on the device or charger. **Use with improper voltage will damage product and VOID warranty.**
- If product is opened/disassembled for ANY reason, warranty is VOID.
- Do not use this product if it has a damaged cord or plug, if it is not working properly, if it has been dropped or damaged, or dropped into water. Contact Waterpik International, Inc. or your local distributor which can be found at www.waterpik.com/intl.
- Do not direct water under the tongue, into the ear, nose or other delicate areas. This product is capable of producing pressures that may cause serious damage in these areas.
- Use this product only as indicated in these instructions or as recommended by your dental professional.
- Only use tips and accessories that are recommended by Water Pik, Inc.

- Do not drop or insert any foreign object into any opening on the product.
- Keep charging unit and cord away from heated surfaces.
- Do not use this product outdoors or operate where aerosol (spray) products are being used or where oxygen is being administered.
- Do not use iodine, saline (salt solution), or water insoluble concentrated essential oils in this product. Use of these can reduce performance and will shorten the life of the product. Use of these products may void warranty.
- Fill reservoir with water or other dental professional recommended solution only.
- Remove any oral jewelry prior to use of this product.
- Do not use if you have an open wound on your tongue or in your mouth.
- If your physician advised you to receive antibiotic premedication before dental procedures, you should consult your dentist before using this instrument or any other oral hygiene aid.
- Close supervision is necessary when this product is used by, on or near children or invalids.

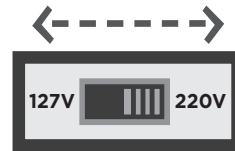- Children should be supervised to ensure that they do not play with the appliance.
- **Your charger can be switched to operate in either 127V or 220V.**
- **The charger is shipped from the factory with the switch in the 220V position.**
- **Before plugging the charger into the wall receptacle, check to ensure the switch is in the correct voltage position for your location.**

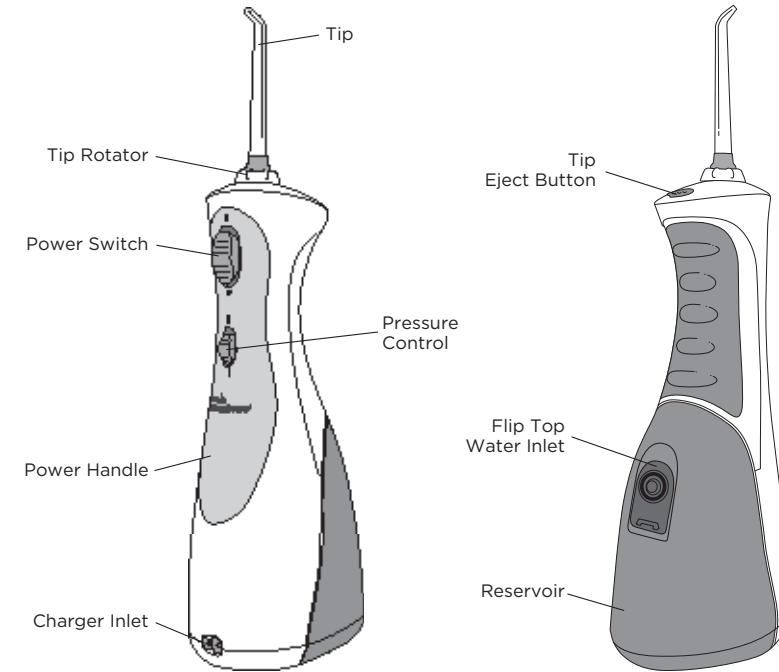
SAVE THESE INSTRUCTIONS.

NOTE: All tips not included in all models.

Replace Every 3 Months	
Orthodontic Tip	BEST FOR
	<input checked="" type="checkbox"/> Braces <input checked="" type="checkbox"/> General Use
Pik Pocket™ Tip	BEST FOR
	<input checked="" type="checkbox"/> Periodontal Pockets <input checked="" type="checkbox"/> Furcations
Plaque Seeker® Tip	BEST FOR
	<input checked="" type="checkbox"/> Implants <input checked="" type="checkbox"/> Bridges <input checked="" type="checkbox"/> Crowns <input checked="" type="checkbox"/> General Use <input checked="" type="checkbox"/> Retainers
Toothbrush Tip	BEST FOR
	<input checked="" type="checkbox"/> General Use

Replace Every 6 Months	
Classic Jet Tip	BEST FOR
	<input checked="" type="checkbox"/> General Use
Tongue Cleaner	BEST FOR
	<input checked="" type="checkbox"/> Fresher Breath

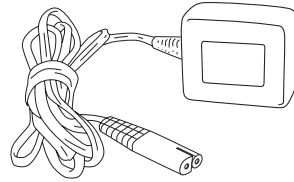Charger
 Waterpik® p/n 20016190 (White)
 Waterpik® p/n xxxxxxxx (Black)

Replacement tips/attachments may be purchased from your local distributor, which can be found at www.waterpik.com/intl.

Charging Unit Prior to First Use

Plug the cord into any standard electrical outlet. Make sure you push the switch into the "OFF" (0) position before charging. The unit will not charge if the switch is in the "ON" (I) position. See "Cleaning and Troubleshooting Your Waterpik® Water Flosser" for guidance on how to look after your battery. The Waterpik® Water Flosser will stay charged for up to 2 weeks depending on use.


Filling the Reservoir

Disconnect charging cord from the unit prior to filling the reservoir and product use. Lift flip top on reservoir and fill the reservoir with lukewarm water.

Inserting and Removing Tips


Insert tip into the center of the knob at the top of the Water Flosser handle. Press firmly until the tip clicks into place. To remove the tip from the handle, press the tip eject button and pull tip from the handle.

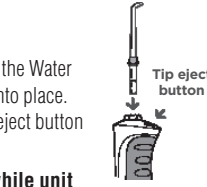CAUTION: Do not press tip eject button while unit is running.

Adjusting the Pressure Setting

You can adjust the pressure setting on your Water Flosser by simply sliding the pressure switch on the hand

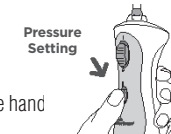Recommended Technique

Lean low over sink and place the tip in mouth. Aim the tip toward teeth, and turn the unit ON (I).


For best results, start with the back teeth and work toward the front teeth. Glide tip along gumline and pause briefly between teeth. Continue until you have cleaned the inside and outside of both the upper and lower teeth.

Direct the jet stream at a 90-degree angle to your gumline. **Slightly close lips to avoid splashing** but allow water to flow freely from mouth into the sink. Keep unit upright during use for best results.

When Finished


Turn the unit OFF (0). Empty any liquid left in the reservoir. If desired, the reservoir may be removed from the power handle by pulling the reservoir release latch while sliding the reservoir down towards the base of the unit. To replace the reservoir on the power handle, simply slide it up towards the top of the unit. The reservoir will snap into place when fully engaged.

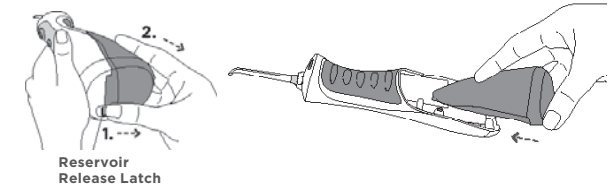TIP USAGE**Pik Pocket™ Tip**

The Pik Pocket™ tip is specifically designed to deliver water or anti-bacterial solutions deep into periodontal pockets.


To use the Pik Pocket™ tip, **set the unit to the lowest pressure setting**. Place the soft tip against a tooth at a 45-degree angle and gently place the tip under the gumline, into the pocket. Turn the unit on and continue tracing along the gumline.

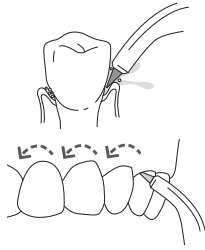**Tongue Cleaner Tip**

To use the tongue cleaner tip, set to lowest pressure setting, place tip in the center/middle of your tongue about half way back and start the Water Flosser. Pull forward with light pressure. Increase pressure as you prefer.

**Orthodontic Tip**

To use the orthodontic tip, gently glide tip along gumline, pausing briefly to lightly brush area between teeth and all around orthodontic bracket, before proceeding to the next tooth.

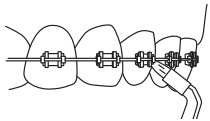**Toothbrush Tip**

Place the toothbrush tip in mouth with the brush head on the tooth at the gumline.

The toothbrush tip can be used with or without toothpaste. Turn the Water Flosser on so the water flows through the tip. Using a **light pressure** (bristles should not bend) massage the brush back and forth with very short strokes-much like you would with a manual toothbrush.

Plaque Seeker® Tip

To use the Plaque Seeker® tip, place the tip close to the teeth so that the bristles are gently touching the teeth. Gently glide tip along gumline, pausing briefly between teeth to gently brush and allow the water to flow between the teeth.

**Using Mouthwash and Other Solutions**

Your Waterpik® Water Flosser can be used to deliver mouthwash and antibacterial solutions. After using any special solution, rinse unit to prevent clogging by partially filling reservoir with warm water and running unit with tip pointed into sink until unit is empty.

Cleaning

Never immerse unit in water. Before cleaning, unplug from the electrical outlet. Clean the product when necessary by using a soft cloth and mild non-abrasive cleanser to wipe exterior.

The water reservoir is removable for easy cleaning and is top rack dishwasher safe. It is recommended that the reservoir is removed and cleaned weekly.

Removing Hard Water Deposits

Hard water deposits may build up in your unit, depending on the mineral content of you water. If left unattended, it can hinder performance. Cleaning internal parts: add 2 tablespoons of white vinegar to a full reservoir of warm water. Point the handle and tip into sink. Turn unit ON and run until reservoir is empty. Rinse by repeating with a full reservoir of clean warm water.

This process should be done every 1 to 3 months to ensure optimal performance.

Battery Care and Disposal

The unit will not charge if the switch is in the ON (I) position. If you use your Water Flosser once a day or less, it is recommended you charge it overnight, once a week. If you use it twice a day or more, charge the unit more frequently. Do not run the batteries down completely as this may shorten the lifespan of the batteries and will require charging the unit for 24 hours to obtain a full charge.

If the Water Flosser is to be stored for an extended period of time (six months or more), be sure to fully charge it prior to storage. The unit contains a Nickel Metal Hydride battery. At the end of the battery life, please recycle the unit at your nearest recycle center per local or state requirements. The battery is not replaceable.

Service Maintenance

Waterpik® Water Flossers have no consumer-serviceable electrical items and do not require routine service maintenance. For all your service needs

go to www.waterpik.com/intl to find your local distributor. Please do not ship product to Waterpik International, Inc. offices, this will delay service.

Some replacement parts and accessories for the Water Flosser are also available from your local distributor which can be found at www.waterpik.com/intl.

Refer to the serial and model numbers in all correspondence. These numbers are located on the back of the handle and bottom of the reservoir respectively.

Still have questions? Go to www.waterpik.com/intl to find your local distributor.

Water Pik, Inc. warrants to the original consumer of this new product that it is free from defects in materials and workmanship for 1 year from the date of purchase. Consumer will be required to submit the original purchase receipt as proof of purchase date and if requested, the entire product, to support a warranty claim. Water Pik, Inc. will replace at its discretion any part of the product, which in its opinion is defective, provided the product has not been abused, misused, altered or damaged after purchase, was used according to instructions, and was used only with accessories or consumable parts approved by Water Pik, Inc. Installation is the responsibility of the consumer and is not covered by the warranty. This limited warranty excludes accessories or consumable parts such as tips, etc. **Use with improper voltage will damage product and VOID warranty.**

This warranty gives you specific legal rights. You may also have other rights that vary, depending on the law in your place of residence.

Water Pik, Inc.

1730 East Prospect Road
Fort Collins, CO 80553-0001 USA

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Form No. 20021794-F AA
©2016 Water Pik, Inc.

WaterPik, International, Inc.

Castle Court, 41 London Road
Reigate, Surrey RH2 9RJ, United Kingdom

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Waterpik® is a trademark of Water Pik, Inc. registered in Argentina, Australia, Benelux, Brazil, Canada, Chile, China, Czech Republic, Denmark, EU, Finland, France, Germany, Hong Kong, Hungary, India, Israel, Italy, Japan, Korea, Mexico, New Zealand, Norway, Pakistan, Poland, Russian Federation, Saudi Arabia, Singapore, South Africa, Spain, Sweden, Switzerland, Taiwan, Turkey, Ukraine, United Kingdom, and the United States.

Waterpik® (stylized) is a trademark of Water Pik, Inc. registered in Australia, Brazil, Canada, China, EU, Hong Kong, Israel, Japan, Korea, Mexico, Russian Federation, Saudi Arabia, South Africa, Switzerland, Taiwan, Trinidad & Tobago, and the United States.

Pik Pocket™ is a trademark of Water Pik, Inc.

Plaque Seeker® is a trademark of Water Pik, Inc. registered in Canada, China, EU, and the United States.

Designed in USA. Made in China.

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Pik Pocket™ é uma marca registrada da Water Pik, Inc.
Plaque Seeker® é uma marca registrada da Water Pik, Inc. no Canadá, China, UE e nos Estados Unidos.
Pik Pocket™ es una marca de comercio de Water Pik, Inc.